



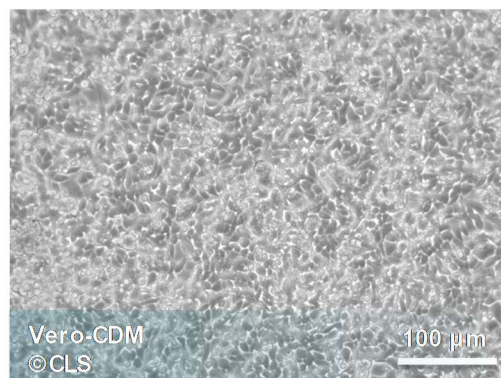

# SciN4Vero

SciN4Vero is a chemically defined, serum- and animal component-free medium supplemented with recombinant human growth factors. The unique formulation is designed to propagate the proprietary Vero-CDM (Art. No. 605372-CDM). The low protein content (<2 mg/L) facilitates recombinant protein and/or virus production.

Adherent Vero-CDM cultured and adapted in SciN4Vero:

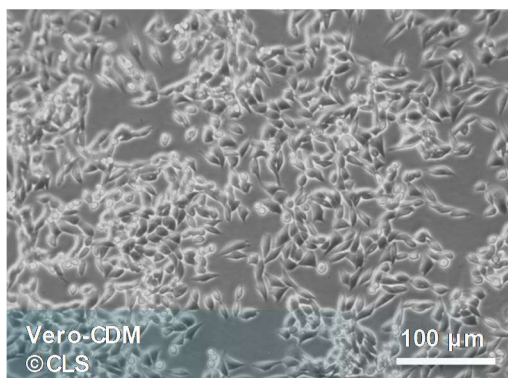

## Ordering Information

Cat-No.	Designation	Size
820502	SciN4Vero Chemically defined, serum-free, animal component-free medium. With L-glutamine.	500 mL
820503	SciN4Vero Chemically defined, serum-free, animal component-free medium. Without L-glutamine.	500 mL

pH = 7.3 +/- 0.2 at 20-25°C.

Each lot of SciN4Vero has been tested for sterility and absence of Mycoplasma and Bacteria.

Stable for about one (1) year from date of manufacture as indicated on the label.

Manufacturer is Scinora GmbH, Heidelberg, Germany. CLS Cell Lines Service GmbH holds the exclusive distribution. Quality control and documentation are performed according to ISO9001:2015 based on GLP.



## Maintenance

- Keep refrigerated at +2°C to +8°C in the dark. Freezing as well as warming up to +37°C minimize the quality of the product.
- Do not heat the medium to more than 37°C or use uncontrollable sources of heat (e.g. microwave appliances).
- If only a part of the medium is to be used, remove this amount from the bottle and warm it up inside a separate device.

**This product is for *in vitro* use and manufacturing only.**

**Not intended for clinical or diagnostic applications.**